LAVA heat pack
Saco Térmico de LAVA
Sac Chauffant LAVA

LAVAcare
PHYSIO≡THERM®

Made in Germany by Physiotherm GmbH
Lavacare - Registered Trademark

Instructions for use

Designed in consultation with physicians, chiropractors and physical therapists

Lavacare is not an FDA-approved device, and intended solely for home therapy.
Consult your physician if you have any questions as to the medical or health benefits of this device.

Operating principle

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Principe d'action

EN Heat is an ancient remedy, which LAVAcare puts to good use. The pack is filled with special lava sand, which is characterized by its excellent ability to store heat. The lava mineral particles inside the pack give off a deep heat that is particularly pleasant, effective and completely dry. The heat also could reach deeper-lying tissue, may relieve pain and relax muscles, may promote circulation and may stimulate metabolism. LAVAcare delivers long-lasting, comforting warmth and adapts to the body's own heat after use. In everyday use, LAVAcare is also used for preheating beds, wheelchairs and car seats. LAVAcare is also used in physical therapy as an aid to improve the movement of individual joints with heat, massage and special exercise.

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
Note: As a condition of sale, the purchaser assumes responsibility for proper care and use of this product in accordance with our printed instructions. The purchaser or user must be the judge of when to use the product and for how long.

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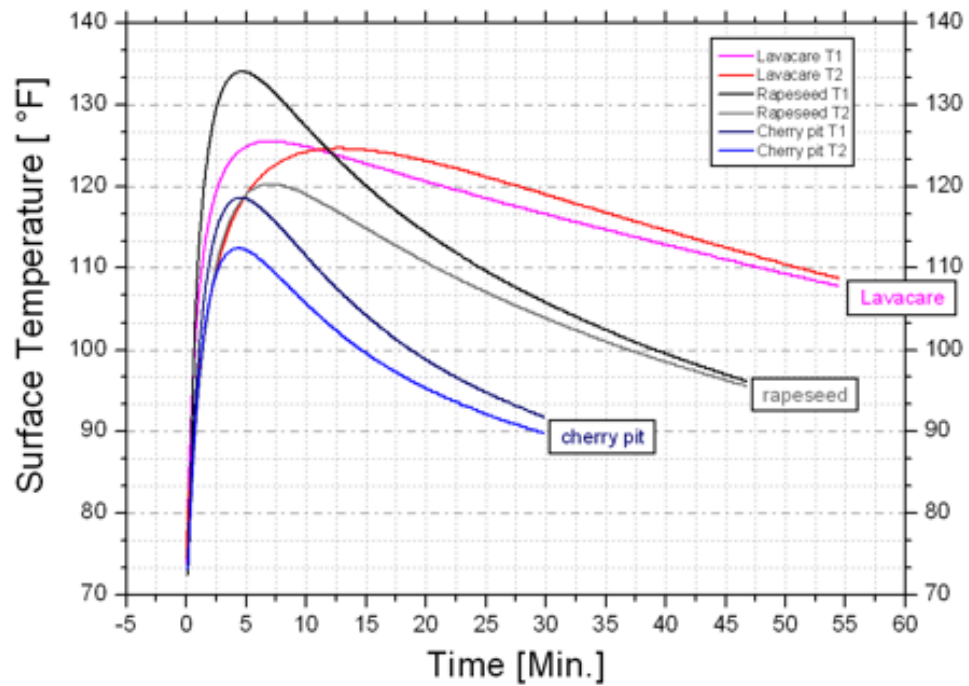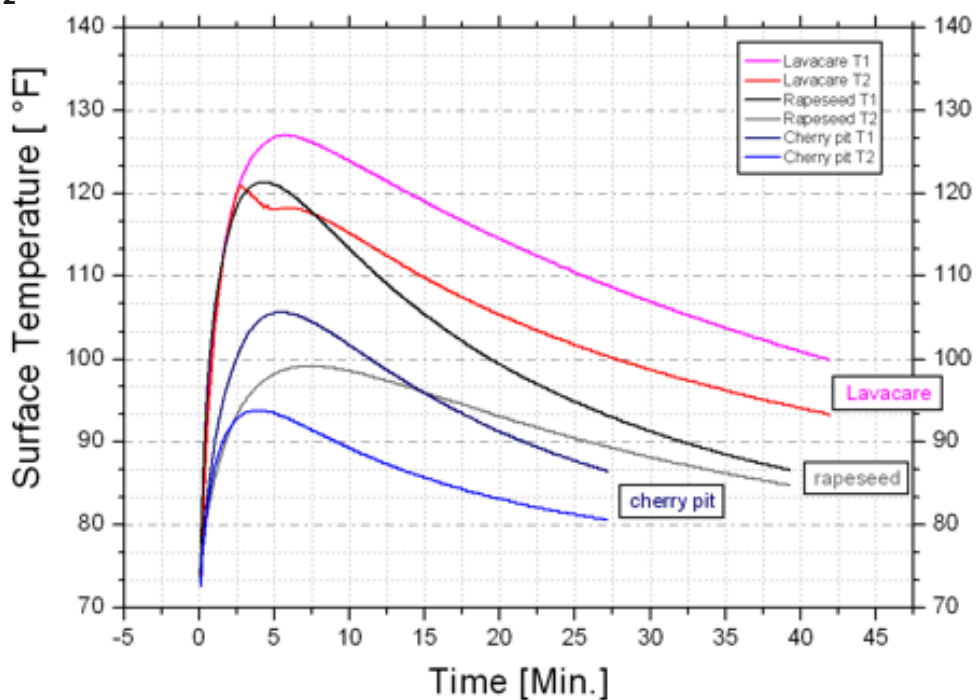
EN Misuse or abuse may damage or destroy this product. Read and follow all instructions carefully. LAVACARE is guaranteed for a period of 1 year from the date of purchase against failure of the inner pack or the patented LAVACARE compound. Warranty is void if unit has not been used according to printed instructions, or if pack is broken or shows evidence of overheating. LAVACARE will be replaced at no cost during the warranty period.

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Note: Valid only in the USA

Thank you for purchasing LAVACARE products
Please fill this card and return it to validate your warranty

Name _____

Address _____

City _____ State _____ Zip _____

Age _____ ☐ Male ☐ Female

Purchase price \$ _____ Product size: ☐ 12" x 8.25"

Where was this purchase made? _____

How did you hear about LAVACARE? _____

On what area of your body do you use your LAVACARE? _____

Comments _____

All information is confidential. Ditum Group LLC does not sell information to other business

Important Information Información Importante Informations Importantes

EN Caution:

- Follow instructions carefully. Use LAVA CARE only as personal therapeutic pack. Misuse can cause burns. Always use the pack cover.
- Test the temperature before using, by wrapping pack in a towel, and applying LAVA CARE to sensitive skin.
- Do not overheat. Overheating in the microwave can damage or destroy the product.
- If the product starts to inflate while heating, immediately turn off the microwave and wait for the inflation to go down before opening the door.
- Due to a higher power of the microwave oven or a longer heating time than indicated in the instructions for use the fabric develops brown stains. The brown stains are caused by the flame retardant treatment of the pack. Overheating of the pack will create black stains and ruptures which can produce smoke and leakage of the lava sand.
DO NOT OVERHEAT THE PRODUCT – IT WILL BE DAMAGED. Damaged packs can not be reheated or used again.
- Persons with diabetes, circulatory problems, nerve damage, paralysis or sensitive skin should use only as directed by a physician. Consult your doctor before treating any serious injury or if pain persists.
- Do not use the LAVA CARE in conjunction with liniments, balms or lotions.

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For more information visit our website: www.lavacare.us and www.lavacare.ca

Place
Postage
Here

Ditum Group, LLC

205 Pecan Blvd
McAllen, TX 78501-2354
United States
Ph. 1-800-392-6675
Fax: 1-866-237-0713

Important Information Información Importante Informations Importantes

- EN** The heat from the LAVA CARE heat pack is also used to help:
- Relieve tense muscles, back, joint, stomach and abdominal pain.
 - The heat may reach deep tissue areas, may alleviate pains.
 - Could improve local circulation (e.g. cold feet).

**Note: Always consult your doctor if unsure about the cause of your pain.*

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Lavacare is not an FDA-approved device, and intended solely for home therapy.
Consult your physician if you have any questions as to the medical or health benefits of this device.

USA and Canada Distributor: Ditung Group, LLC - 205 Pecan Blvd, McAllen, TX 78501-2354, United States
Ph. 1-800-392-6675 Fax: 1-866-237-0713 /www.lavacare.us / www.lavacare.ca

Distribuidor en México: Metrofarma S.A de C.V - Isidoro Sepúlveda 780-C, Parque Industrial Multipark,
Apodaca, N.L, CP. 66600, México - Tel. +52 (811) 520-1123/25 - Lada sin costo: 01-800-836-5491 /www.lavacare.com.mx

Instructions for use

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Mode d'emploi

EN Heat the pack in the microwave for **TWO MINUTES (600 WATTS)**
Please ensure that the turntable can turn freely and do not leave the appliance unattended while the pack is being heated!

Warning: Remove the pack cover before heating.

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EN Alternatively heat the pack in an oven at approx. 150° C, 302° F for 15 - 20 minutes.

Warning: Remove pack cover before heating.

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
EN **Caution:** Risk of overheating and fire due to prolonged heating of the pack! To avoid overheating, do not reheat the pack for at least 40 minutes after the initial application.

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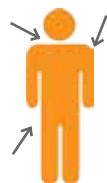EN After heating the pack, place it on the affected part of the body. If your skin is sensitive to temperature, wrap the pack in a towel.

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EN **LAVACARE** can also be used as a cold compress. Place the pack in the freezer compartment for 15 - 20 minutes before use.

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